



sonicaid

Team 3 NST

ADVANCED FETAL MONITORING FOR THE NEXT GENERATION

- Superior Performance for High BMI Moms
- Easy Touchscreen Operation
- Outstanding Performance

Sonicaid Team 3 NST

Fetal monitoring for the next generation

Team 3 NST improves fetal heart rate monitoring with a large, color touchscreen that displays the fetal heart rate in two modes: big numbers and trace view.

Large, high contrast number mode allows the heart rate to be easily seen from a distance. The trace mode allows for review and scrolling on screen for in-depth review.

Along with high sensitivity transducers, ideal for high BMI moms, with unique digital processing, Team 3 NST delivers the best possible fetal monitoring for the next generation.



8.4" COLOR TOUCHSCREEN

- High resolution display can be viewed from a distance even in brightly or dimly lit environments
- Wide viewing angle
- Can be used with gloves on




MULTI FUNCTIONAL DISPLAY

- Toggle between "BIG Numbers" and FHR trace display modes
- On screen trace review and scrolling with clear color separation of multiple traces



OPTIONAL MATERNAL VITAL SIGNS

- NiBP - clinically validated on pregnant moms
- Measures on inflation, increasing comfort and reducing time

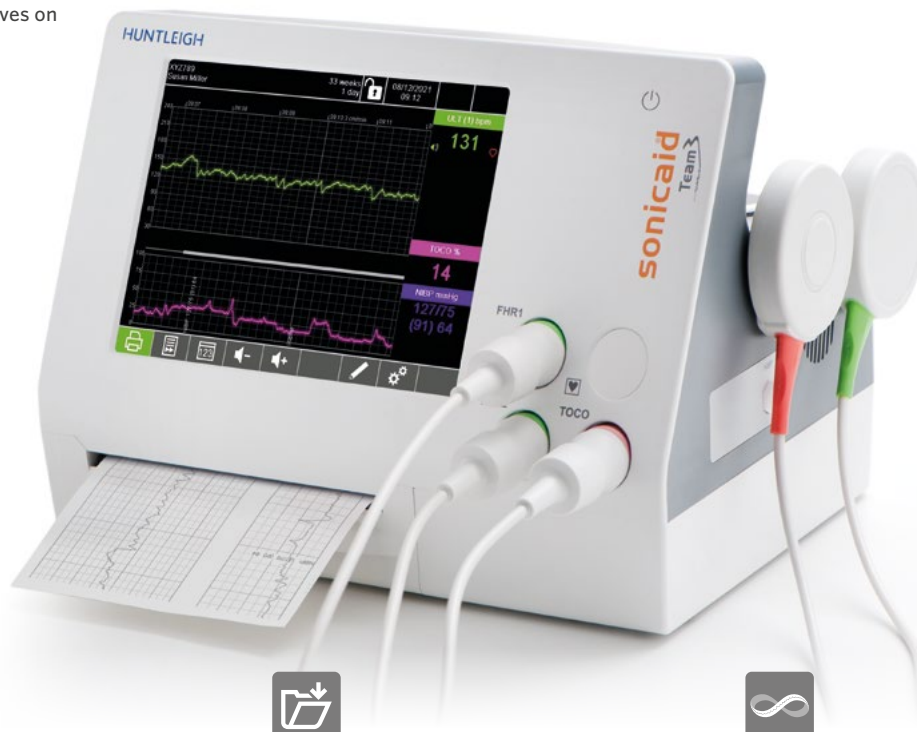

OPTIONAL INTEGRATED BATTERY

- Allows portability for up to 4 hours



PATIENT RECORDS

- Patient database
- On-screen gestational age calculator (LMP, GA, EDD)
- Records stored on a USB drive for easy transfer to the computer



TRACEBANK

- Never lose a trace again by storing fetal heart traces in the monitor's internal memory and backup feature



High BMI, high performance transducers

Widebeam ultrasound technology

Widebeam ultrasound technology provides a wider, deeper, uniform beam shape to quickly and easily locate the heartbeat and maintain contact.

Patented "Locate and Track" technology

Maintains contact with the fetal heart, especially on moms with high BMI, reduces loss of contact and signal.

Managing multiple pregnancies

Team 3 is capable of monitoring single and twins, allowing just one monitor to be used.



Software solutions – archive, review or share

Sonicaid Reporter Software

Sonicaid Reporter allows for reviewing, storing and printing of patient records. Sonicaid Reporter included at no charge.

Downloadable to USB

Fetal heart rate traces can be quickly downloaded on a USB drive and transferred to any computer.

Easy integration into EMR

Traces can be saved as a PDF for simple integration into EMR.



Models								
	TEAM3A-USA	TEAM3A-E-USA	TEAM3A-B-USA	TEAM3A-EB-USA	TEAM3A-N-USA	TEAM3A-EN-USA	TEAM3A-NB/USA	TEAM3A-ENB-USA
8.4" Color Touchscreen	•	•	•	•	•	•	•	•
Internal Trace Memory	•	•	•	•	•	•	•	•
Twins Capable*	•	•	•	•	•	•	•	•
Ultrasound and TOCO Transducer	•	•	•	•	•	•	•	•
Patient Event Marker	•	•	•	•	•	•	•	•
Maternal Vital Signs (NIBP and Heart Rate)					•	•	•	•
Printer	•		•		•		•	
Rechargeable Battery			•	•			•	•
Weight	13.5 lbs (6Kg) Max							
Dimensions (LxWxH)	9"x12.6"x9.4"	9"x12.6"x7.5"	9"x12.6"x9.4"	9"x12.6"x7.5"	9"x12.6"x9.4"	9"x12.6"x7.5"	9"x12.6"x9.4"	9"x12.6"x9.4"

*Twins Capable - supplied with only 1 ultrasound transducer as standard

Technical Specifications	
Interfaces	RS232, USB (x2), Wireless System port
Mains power	85-264V, 50/60Hz, <100VA, IEC power cable

Accessories			
ACC-OBS-008	Team3 / FM800 Encore 1.0MHz Ultrasound Transducer (green) - 6.6 ft	ACC-MI1160	Hose for Blood Pressure Cuff
ACC-OBS-008L	Team3 / FM800 Encore 1.0MHz Ultrasound Transducer (green) - 9.8 ft	SP7775-6901	Event Marker
ACC-OBS-009	Team3 / FM800 Encore Toco Transducer (Pink) - 6.6 ft	ACC3	Ultrasonic Gel (box 12x250ml Bottles)
ACC-OBS-009L	Team3 / FM800 Encore Toco Transducer (Pink) - 9.8 ft	ACC-MI1136	Latex free - Reusable transducer belts (pair)
ACC-OBS-073	Medium Blood Pressure Cuff, Arm Circumference Range 24-32cm	ACC220	Latex free - Disposable Single use transducer belts (pack of 100 pairs)
ACC-OBS-074	Large Blood Pressure Cuff, Arm Circumference Range 32-42cm	ACC-8400-8003	Standard grade chart paper (Box 20 pack)
ACC-OBS-075	Extra-Large Blood Pressure Cuff, Arm Circumference Range 38-46cm	ACC-321414	ArchiTrace 25 year archival chart paper (Pack of 20)

Warranty	
2 years on fetal monitor	
1 year on transducers and accessories	

As a proud member of the Arjo family, we have been committed to supporting healthcare professionals in improving outcomes and enhancing patient wellbeing since 1979. We do this through our proven solutions for Vascular Assessment & Treatment and Fetal & Patient Monitoring. With innovation and customer satisfaction as our guiding principles, we strive for clinical excellence and improved performance, for life.

Distributed in the USA by:

Arjo Inc..
 2349 West Lake Street, Suite 250 Addison, Illinois 60101
 T: +1 (800) 323 1245 F: +1 888 594 2756
 www.huntleigh-healthcare.us

Registered No: 942245 England & Wales, Registered Office:
 ArjoHuntleigh House, Houghton Hall Business Park, Houghton Regis, Bedfordshire, LU5 5XF
 ©Huntleigh Healthcare Limited 2015-2020

A Member of the Arjo Family

As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.

- All models supplied complete with gel, latex free transducer belts, chart paper, user manual and quick start guide.

AW:1001047-4